

NEOCONNECT[®] for the Pharmacy

Where ENFit[®] and Innovation Meet



NEOMED
www.neomedinc.com

NEOCONNECT
with ENFit[®] Connectors



(770) 516-2225
100 Londonderry Ct, Suite 112
Woodstock, GA 30188
www.neomedinc.com

ENFit® transition is underway. Has your Pharmacy identified a supplier that delivers a solution preserving current feeding/medication delivery protocols?

- Can we use the same syringe for oral and enteral delivery?
- Does our pharmacy supplier offer syringe sizes of 5 mL or smaller with the ENFit Low Dose Tip syringe design that complies with the July 2016 GEDSA Guidance Document?
- Can we deliver accurate low volume doses within current practice?
- Does your pharmacy vendor offer a full complement of ENFit accessories including filling devices with seven sizes of bottle caps and dispensing devices like the Dosemate that meet recommendations from GEDSA and ISMP?
- Do we need a complete portfolio of amber barrel syringes in both bulk and sterile?
- Do we want products that offer minimal changes to current processes?
- Are hands-free self-righting tip caps essential to our aseptic techniques?
- Do we have a single vendor that offers a comprehensive solution for pharmacy and NICU?

NeoMed's NeoConnect is the comprehensive solution.



Pharmacy Syringes

- Metric scale gradients only (Complies with ISMP Best Practices #5)¹
- Available in sizes 0.5 mL to 60 mL in clear with purple gradient markings or amber with white gradient markings
- FDA 510(k) clearance allows both enteral and oral administration in all sizes, eliminating the need to stock multiple syringe types²
- Hands-free, self-righting, **click**-to-close tip caps support aseptic technique
- Sterile amber syringes available in sizes 1 to 12 mL

Low Dose Syringe

Designed to minimize process changes in both pharmacy and bedside filling and dispensing operations

- Internal tip lumen mirrors traditional male syringe designs reducing deadspace and the potential dosing error associated with backflow
- 0.5 mL, 1 mL, 3 mL, and 6 mL (clear and amber)
- Demonstrated clinical dosing accuracy up to and including dispensing volumes of 5 mL and smaller in accordance with the July 2016 GEDSA Guidance Document
- Low priming volume/minimal dead space
- FDA 510(k) clearance for both oral and enteral delivery²
- Complies with August 2016 GEDSA Position Statement Supporting ENFit Low Dose Tip Syringe



NeoConnect NeoSecure™ Self-Righting Tip Cap

- Self-rights on any hard surface for hands free placement: ideal for aseptic technique
- Secure **click-to-close** feature facilitates hands-free operation and eliminates the need for twisting
- Cap is secure during transport and refrigeration while helping mitigate the risk of incidental leakage
- Minimizes impact to process flow



DoseMate Oral Administration Tip


- Available in two sizes
Regular: Nipple style device allows for improved patient acceptance during direct oral administration.
Long: The Directional Long (DL) is designed to allow for dosing of meds deeper into the patient's mouth when appropriate
- Ergonomically designed for patient comfort during direct oral administration without introducing a choking hazard⁶
- Addresses concerns with threaded syringe tips or safety tabs during oral administration
- Both provide an easy device solution that follows the August 2016 GEDSA Guidance Document recommendation to clear out the LDT syringe "moat" (the area between the male lumen and the outer ring that lies outside of the fluid path)

NeoMed has over 100 active and pending patents protecting our product portfolio.



Pharmacy Caps

- Promotes dosing accuracy by displacing the same volume as the ENFit male connector
- Sizes 18 mm to 38 mm
- Designed for rapid filling operations
- Non-threaded design helps minimize impact to process flow
- Tethered cap for secure closure

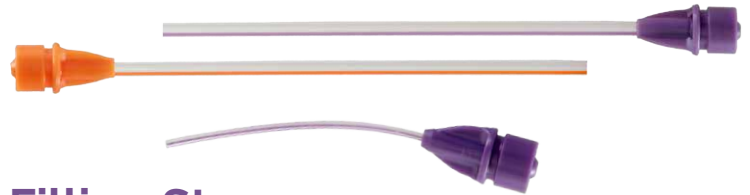

Filling Straws

- Promotes aseptic technique
- Can be used at bedside
- Small bore and large bore available



Syringe to Syringe Coupler

- Connects two NeoConnect syringes for transfer or mixing
- Facilitates patient-specific dosing from another NeoConnect syringe

¹ Institute for Safe Medication Practices. 2014–15 Targeted medication safety best practices for hospitals. (Available at: www.ismp.org/Tools/BestPractices/default.asp).

² Food and Drug Administration Department of Health & Human Services. "NeoConnect Oral/Enteral Syringes with ENFit Connector (12 mL to 100 mL) And NeoConnect Low Dose Tip Oral/Enteral Syringes with ENFit Connector (0.5 mL to 6 mL)." K161039. June 2016.

³ Method For Identifying Toys And Other Articles Intended For Use By Children Under 3 Years Of Age Which Present Choking, Aspiration, Or Ingestion Hazards Because Of Small Parts.

Ordering Information

Sterile Individually Packaged NeoConnect Pharmacy Syringes Amber Barrel with White Gradient Markings			
Product Code	Description	UOM	QTY
BA-S1NC-S	1 mL	CA	200
BA-S3NC-S	3 mL	CA	200
BA-S6NC-S	6 mL	CA	200
BA-S12NC-S	12 mL	CA	200

Bulk Packaged NeoConnect Pharmacy Syringes Amber Barrel with White Gradient Markings			
Product Code	Description	UOM	QTY
BA-S05NC	0.5 mL	CA	500
BA-S1NC	1 mL	CA	500
BA-S3NC	3 mL	CA	500
BA-S6NC	6 mL	CA	500
BA-S12NC	12 mL	CA	500
BA-S20NC	20 mL	CA	200
BA-S35NC	35 mL	CA	200
BA-S60NC	60 mL	CA	200

Pharmacy Bottle Caps			
Product Codes	Description	UOM	QTY
PC-ANC	Size A, 18 mm	CA	25
PC-BNC	Size B, 20 mm	CA	25
PC-CNC	Size C, 22 mm	CA	25
PC-DNC	Size D, 24 mm	CA	25
PC-ENC	Size E, 28 mm	CA	25
PC-GNC	Size G, 33 mm	CA	25
PC-INC	Size I, 38 mm	CA	25

NeoConnect Pharmacy Accessories			
Description	Product Code	UOM	QTY
NeoConnect™ DoseMate™ Oral Administration Tip	POA-TNC	CA	50
NeoConnect™ DoseMate™ DL (Directional Long) Oral Administration Tip	POA-TNCL	CA	50
Syringe-to-Syringe Coupler, ENFit™ Syringe Compatible, Sterile	PCO-FNC	CA	100
Syringe-to-Syringe Coupler, ENFit™ Syringe Compatible, Non-Sterile	PCO-FNC-NS	CA	400
NeoConnect™ NeoSecure™ Tip Cap, Sterile	PNM-TCNC	CA	200
NeoConnect™ NeoSecure™ Tip Cap, Non-Sterile	PNS-TCNC	CA	600
Transfer Lid, NeoConnect™ Purple	PTL-NC	CA	200

Bulk Packaged NeoConnect Pharmacy Syringes Clear Barrel with Purple Gradient Markings			
Product Code	Description	UOM	QTY
BP-S05NC	0.5 mL	CA	500
BP-S1NC	1 mL	CA	500
BP-S3NC	3 mL	CA	500
BP-S6NC	6 mL	CA	500
BP-S12NC	12 mL	CA	500
BP-S20NC	20 mL	CA	200
BP-S35NC	35 mL	CA	200
BP-S60NC	60 mL	CA	200

NeoConnect Straws with ENFit Connector (50/case)		
Orange	Purple	Description
EXT-3STNC	PEXT-3STNC	3" (8 cm) PVC Straw
EXT-6STBNC	PEXT-6STBNC	6" (15 cm) PVC Straw

All 0.5, 1, 3, and 6 mL syringes have the low dose tip syringe design.